



# Amfibi™ Basic

Hygiene Commode



# Amfibi™ Basic

## Hygiene Commode

**Basic** was developed to meet most basic needs. With its attractive price it provides a lot of hygiene chair at low cost. **Basic** has adjustable seat height and a wide range of accessories that make it easy to individually adapt the chair to suit every patient's needs.

### Flip-up armrests

Armrests are easy to flip-up.

### Different seat alternatives

There are several different seat alternatives available. Use the soft PU seat for greatest comfort. See *accessories*.

### Guide included

The guide provided with the chair fits both the Basic bedpan and metal hospital care bed pans.

### Flip-up footrests

The footrests flip up and are easy to swing aside. Can be removed if necessary.

### Durable materials

Manufactured from corrosion-proof materials.

### Comfortable

The backrest is comfortable, adjustable and machine washable at 70 degrees.

### Adjustable height

Telescoping legs easily adjust for seat heights from 49 to 62 cm.

### Lockable castors

All four castors are lockable.




## Technical specifications

Dimensions apply to walkers with 125mm castors

CE marked according to MDD 93/42/EEC

Product Code	Product	Width	Max user weight	Seat Depth	Seat Height	Armrests/internal dimension
RHGA-56-280	Amfibi Basic	57cm	150kg	47cm	49-62cm	47cm
RHGA-56-280XW	Amfibi Basic*	62cm	150kg	52cm	49-62cm	52cm

\* Extra Wide

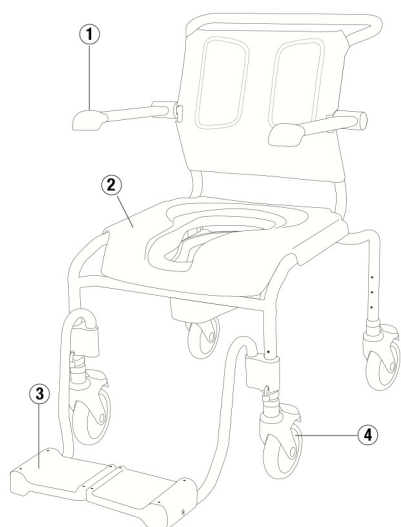

### Maintenance

The commode is made of steel tubes (outer/inner coated) and aluminium. Exposed details that haven't been coated are made of aluminium/stainless steel. The commode can be cleaned with normal cleaning agents (could be cleaned in an autoclave 85 degrees). Check that the chair is not damaged and there are no loose parts.

### Regular checks

1. Armrest- ensure that the glide protection/end caps are free from damage and fitted correctly.
2. Seat/softseat- is complete and without damage. There should be no sharp edges.
3. Footrest- is complete and can be adjusted. Make sure that the clips are complete and tightened.
4. Wheel- should spin easily. Check that the wheels are correctly tightened in the frame.

## Accessories



**Soft Seat**  
RHGA-56-206



**Heel Loops**  
RHGA-56-220



**Safety Bar**  
RHGA-56-212



**Safety Bar for DB**  
RHGA-56-212-DOUBLE



**Side Supports**  
RHGA-56-216




**Splash Guard**  
RHGA-56-213



**Armrest Locking**  
RHGA-56-214



**Waste Container**  
RHGA-56-209

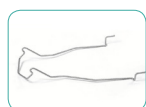

**Waste Container Holder**  
RHGA-56-218



**Safety Belt Hip**  
RHGA-56-217



**Safety Belt Chest**  
RHGA-56-229



**Amputee Support**  
RHGA-56-289-L/R



**Calf Support**  
RHGA-56-288-L/R



**Calf Strap**  
RHGA-56-225



**Castors with brake**  
(various sizes)

SHI200115-v2.2

## Designed to meet needs

GATE hygiene chairs have a unique pillar design for quick and easy seat height adjustment. There is a choice of manual or electric adjustment to suit various needs. Amfibi fits most people no matter how tall, and staff will not need to waste time and energy looking for a suitable hygiene chair for each situation.

## Everyday life improved

And it goes without saying that Amfibi never compromises on user comfort and personnel ergonomics. Individually adaptable comfort and function makes the hygiene chair a much more

pleasant experience for patient and healthcare staff alike. Amfibi is manufactured from corrosion-proof materials and its actuator can withstand high pressure hot water (IP66).

## Best possible comfort

Users often sit in a chair for a long while when washing or taking a shower, and this can become very trying. Amfibi was created to provide extremely high comfort. The seated Amfibi user is able to assume a position for the best possible comfort. The hygiene chair quickly adapts to a personal fit with the aid of easily operated settings, and all models are height-adjustable.

The Amfibi's seat and back are adapted for an upright posture that does not close the pelvis. The ergonomic design minimises pressure on the hips, tendons and blood vessels, while its rounded front edge contributes to a comfortable leg position.

The new, soft, water-repellent backrest that provides the best possible comfort and hygiene is unique to Amfibi (accessory). It is easily adjustable for best support and seat depth.

GATE hygiene chairs are available in three different versions. Amfibi Double is our most advanced, versatile hygiene chair; it is adjustable in height and tilt angle of inclination. Amfibi XL was designed for the bigger, heavier patient, while Basic is our simplest version.

## Other Commodes in the Amfibi™ range



**Amfibi XL** was designed for the bigger, heavier patient and therefore provides extra space and strength and electrical height adjustment. Maximum user weight is 250 kg.

**Amfibi™**  
XL Hygiene Commode



**Amfibi Double** is GATE's most advanced and versatile hygiene chair. It is equipped with twin electric drives, electrically operated raising and lowering and electrically operated tilt. The chair can be adjusted steplessly even when the patient is in the chair.

**Amfibi™**  
Double Hygiene Commode



**Amfibi Double HighBack** is a further development of our much appreciated Amfibi concept. It is equipped with a high, upholstered back with upholstered back including neck rest and integrated seat with rear opening.

**Amfibi™**  
Double HighBack  
Hygiene Commode

## Also Available



**Bure Rise & Go** is an all-round solution that saves space, time and money. Instead of purchasing both hoists and walkers it is probably often sufficient to acquire the flexible, cost effective **Bure Rise & Go**.

**Bure™**  
Rise & Go™



Apart from being an excellent aid for getting patients up and about, the **StandTall Walker** is also a valuable rehabilitation tool. Rehabilitation using the **StandTall Walker** keeps the patient's muscles and alimentary canal active and contributes to improving the general condition of the otherwise bedridden.

**Bure™**  
StandTall Walkers



Distributed by:

For more information, or to find your nearest dealer, call us on 1300 791 404 or visit [www.gaterd.com.au](http://www.gaterd.com.au)

