



**For Home & Business**  
*Self or Assisted Operation*

*Homes*

*Adult Communities*


*Senior Centers*

*Assisted Living*


**The People Picker Upper!**  
**HFL-300/HFL-400**



*Self-Operation*

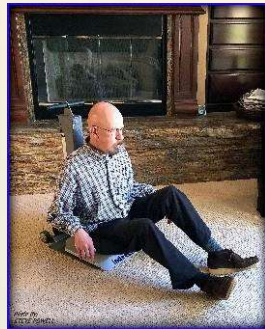

*Assisted Operation*

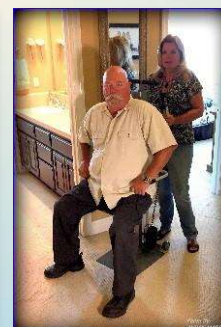



*Self or Assisted Operation!*

*Homes  
Adult Communities  
Senior Centers  
Assisted Living*

## The People Picker Upper! HFL-300/400

### Purpose

The HFL-300 and 400 for Home & Business are designed to lift a fallen person from the floor. After a fall, a person can be easily raised without the risk of injury to themselves or a helper! In the home, the HFL allows for fall recovery without outside assistance, allowing mobility challenged individuals to remain in their own homes longer, thereby maintaining their **privacy, dignity and independence**.

### Operation

Self Operation: After a fall, the fallen person scoots to the HFL, slides onto the beveled seat, pushes a button on the wired remote, activating the motor which raises them from the floor.

Assisted Operation: A helper can easily wheel the HFL to the fallen person and assist them onto the seat, if necessary.



### Physical Specifications

#### **HFL-300 HFL-400**

Safe Working Load .....	<b>300 lbs.</b>	<b>400 lbs.</b>
Maximum Length .....	21"	21"
Maximum Height .....	52"	52"
External Width .....	15"	15"
Maximum Seat Height .....	21"	21"
Wheels .....	5"	5"
Power Pack .....	1.5lbs	2.2 lbs.
Total Weight .....	70 lbs.	72 lbs.

### Summary of Features

- Lifts a person weighing up to 300/400 lbs. to a height of 21" above the floor
- Beveled seat-front for easy slide-on access
- Wired remote for convenient operation
- Swinging/removable Rise Assist Handles for easy transfer to wheelchair
- Self or assisted operation
- Lightweight – 52/56 lbs.
- Small footprint for tight spaces
- Powder coated for easy cleaning
- Powered by a rechargeable Lithium Ion battery pack
- Battery pack and charger included

### Electrical Specifications

Standard 2-prong 110V AC Power 9' cord  
Operating Environment: >+5°C <+40°C  
Lithium Ion Battery: 29.9V DC 1800 mHa

*Maintain Privacy, Dignity and Independence*